


Livia

3,5 Fr Urinary Catheter


for the use on premature babies

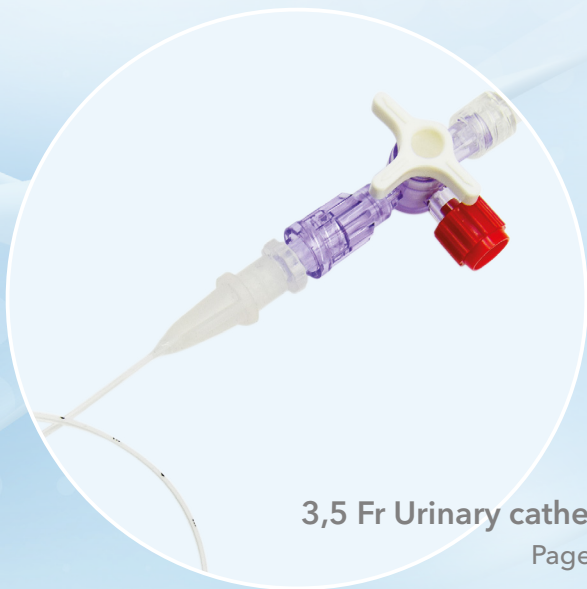

3,5 Fr Urinary catheter
Page 02



Catheter set
Page 03

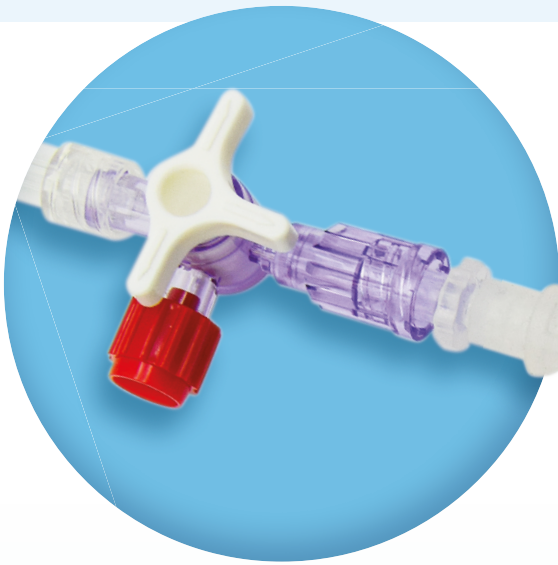

Livia 3,5 Fr Urinary Catheter




The 3,5 Fr Livia urinary catheter (Ø1,2 mm) is a specially developed medical product designed for the use in premature babies. The urinary catheter is used to precisely measure the premature baby's urine output and ensure continuous monitoring.

The 3,5 Fr Livia urinary catheter is used to obtain urine cultures and accurate urine output in neonates. The silicone material of the Livia urinary catheter minimizes urethral irritation and incrustations. This preferred urinary drainage method eliminates the inconvenience and measurement errors associated with weighing diapers.

The 3,5 Fr urinary catheter is a closed system that is used for precise urine excretion in premature babies.

Features:

The silicone catheter has a soft tip to reduce trauma during insertion.

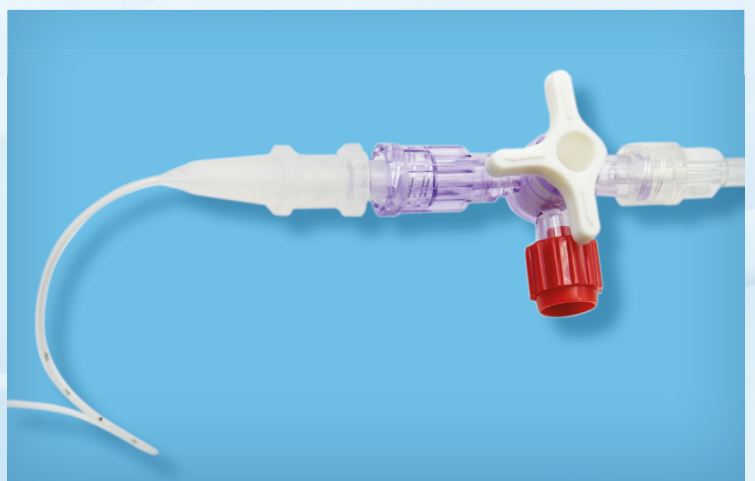

Art. no.	Catheter Ø	Description	Pcs./Pkg.
LK35Fr250	3,5 Fr 1,2 mm	Disposable product, sterile	10

Livia 3,5 Fr Urinary Catheter Set



The Livia catheter set is a closed urinary system that can be used to accurately measure urine out-put in premature babies. The Livia catheter set consists of a preassembled 3.5 Fr silicone urinary catheter, a three-way stopcock for sampling purposes and a 100 ml burette for accurate measurement of urine output.

The pre-assembled, sterile, closed system reduces the risk of infection. The stopcock with injection point enables the sampling of uncontaminated samples in the middle stream. The 100 ml collection burette enables precise measurement of urine excretion. The silicone bladder catheter has a soft tip to reduce trauma during insertion. The catheter is also available separately.

The set consists of:

- biocompatible 3.5 Fr urinary catheter.
- soft, rounded tip.
- stopcock with latex-free injection point.
- 100 ml collection burette with drain tap.
- 100 ml urine collection bag.

Art. no.	Catheter Ø	Description	Pcs./Pkg.
LK35Fr250-Set	3,5 Fr 1,2 mm	Disposable product, sterile	10